




The
precision
expert


THERMOPLASTIC ACTIVATION

Thickness: 1.2 mm, 1.6 mm, 2.0 mm and 2.4 masks


Thickness: 3.2 mm and hybrid masks




Activation temperature:
between 65°C and 70°C
(between 149°F and 158°F)


Activation temperature:
between 65°C and 70°C
(between 149°F and 158°F)




Heating time
min. 3 minutes
max. 30 minutes


Heating time
min. 4 minutes
max. 30 minutes




Drying:
max. 10 seconds


Drying:
max. 10 seconds



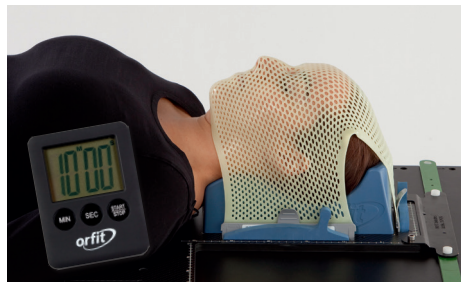

Modelling time:
1 to 1.5 minutes

Modelling time:
1 to 1.5 minutes



Hardening time:
minimum 10 minutes

Hardening time:
minimum 10 minutes